

### **Specifications**


Product Code		221626		
Мо	del	CF302L-A		
Cir	culation type	Closed circulation		
Cir	culating water	Tap water, nybrine (liquid temperature below 10°C)		
Ор	eration temperature range	5°C to 35°C		
Performance	Temperature control range*1	−20°C to 30°C		
	Temperature adjustment accuracy	JTM K05±1.0°C (≧0°C) JTM K05±1.5°C (<0°C)		
	Temperature fluctuation	JIS 3.0°C (≧0°C) JIS 4.0°C (<0°C)		
	Cooling capacity	Approx. 420W at 10°C liquid temp. Approx. 280W at −10°C liquid temp.		
	Max. flow rate (60Hz)	Approx. 10 L/min		
	Max. lift (60Hz)	Approx. 5.7m		
	Temperature control	Refrigerator ON/OFF control		
	Temperature sensor	Pt100Ω		
ion	Temperature setting/display	White LED digital display, key input, minimum digit 1°C		
ırat	Refrigerator and refrigerant	Air-cooling 450W, R404A		
figu		One-touch connector (Flexible rotary L type) Top rear, 1 system		
Configuration	External circulation nozzle	Outer diameter $\phi$ 10 mm rigid hose connection (Inner diameter $\phi$ 9 mm flexible hose mountable with attached hose nozzle		
	Water circulating pump	Magnet pump 15W (60Hz)		
Safety features		Electric leakage breaker, temperature sensor abnormality, Temperature upper/lower limit alarm, Temperature upper/lower limit abnormality, Refrigerator overload relay, Fan motor protection, Pump thermal protector, Refrigerator protection delay timer		
Other features		Drainage hose, Circulation port condensation water drainage hose, Intake port dustproof filter, Cooling operation switch, Circulation pump switch, Calibration offset, Power failure recovery function (selected)		
Standard	Inner tank material	Stainless steel		
	Inner tank capacity	Approx. 3.9L (liquid volume 3.5L)		
	External dimensions (W×D×Hmm)	205×396×535		
	Weight	Approx. 31 kg		
	Power source	AC115V 7.5A (10A) Single-phase		
Accessories		Two hose nozzles (for flexible hose connection), One condensation waterdrainage hose 1m (inner diameter \$\phi\$12 mm Two clamps for hose fixing		

<sup>\*1</sup> This model is not equipped with the heating function. Depending on the ambient temperature and connection conditions, it may not be able toreach -20°C.

#### Performance Curve

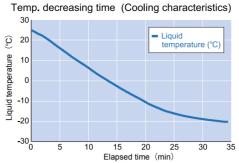

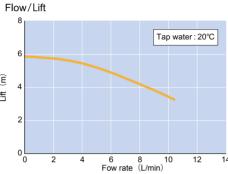

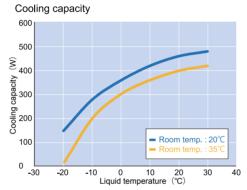

### **Optional Items**

No.	Product Name	Specificatioons	Model	Product Code
1	Heat insulated circulation hose	Inner diameter 9mm, ethylene propylene (soft) hose (2m) 2 pcs., clamps 4pcs.	OCF12	221581
2	Heat insulated circulation hose	High heat insulation, inner diameter 9mm, silicone (soft) hose (2m) 2pcs., clamps 4pcs.	OCF62	221599
3	Heat insulated circulation hose	Outer diameter 10mm, polyurethane (hard) tube (2m) 2pcs.	OCF78	281475
4	L shaped circulation nozzle	φ8mm R3/8 SUS304	OCF80	281476
(5)	L shaped circulation nozzle	φ10.5mm R3/8 SUS304	OCL30	221396
6	L shaped circulation nozzle	φ13mm R3/8 SUS304	OCL40	221397
7	L shaped circulation nozzle	φ16mm R3/8 SUS304	OCL50	221398
8	Flow control valve	R3/8	OCF82	281477
9	External linkage input terminal*2	Operation/stop by no-voltage A contact input	OCF96	281484
10	One touch connector	Straight one-touch connector R3/8×outer diameter \$\phi\$ 10mm for hard tubes	OCF88	281480
(11)	One touch connector	L-shaped swivel type one-touch connector R3/8 × outer diameter $\phi$ 10mm for hard tubes	OCF90	281481
12	Strainer set	Y-type strainer (20 mesh, Rc3/8) +Outlet connection tube	OCF92	281482
13	Hose nozzle	Connection from outer diameter 10mm hard tube to inner diameter 9mm soft hose*3	OCF94	281483

<sup>\*2</sup> Please specify the external linkage input terminal when ordering the main unit. 
\*3 Standard accessories



### For the development of scientific technologies

### **Yamato Scientific Co., Ltd.**

www.yamato-scientific.com

Yamato Scientific Co.,Ltd.

International Sales Department:

Harumi Island Triton Square Office Tower Y, 36F

1-8-11 Harumi, Chuo-ku, Tokyo 104-0053, Japan

TEL: +81-3-5548-7122 FAX: +81-3-5548-0132

<Customer Service/Technical Support> english-website@yamato-net.co.jp

Yamato Scientific Shanghai Corp.

Room 1001-1002, Block B, Xinyan Building, No.65 Guiqing Road, Xuhui District, Shanghai, China TEL: +86-21-6443-5319 FAX: +86-21-5452-0268

TEL: +86-21-6443-5319 FAX: +86-21-5452-0268 <Customer Service/Technical Support> info@yamato-shanghai.com Yamato Scientific America Inc.

925 Walsh Ave.Santa Clara, CA 95050, U.S.A. TEL: +1-408-235-7725 FAX: +1-408-235-7730

<Customer Service Customer> Service@yamato-usa.com

<Technical Support > www.yamato-usa.com





## **Cooling Water Circulator**

(Externally Closed Circulation)
Neocool® Circulator

CF302L-A

# Amazing Price, Performance, Design

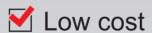

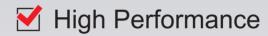

✓ Design • Operability



CF302L-A

Temperature control range: -20°C to 30°C

 $\circ$  Neocool is a registered trademark of Yamato Science Co., Ltd.





### CF302L-A

## **Cooling Water Circulator**

(Externally Closed Circulation)

## CF302L-A

Cooling of the heat generating parts by closed circulation in various applications such as research equipment, industrial machinary, analysis and measuring instruments!







### Easy maintenance of the inlet filter

CF302L-A

Temperature control range: -20°C to 30°C

The front filter mounting plate is only fixed by a magnet, so it can be easily removed by pulling for cleaning. The sheet filter is easy to clean and remount.



### New user friendly controller

A new controller with high visibility and improved operability with a high brightness white LED is adopted. Switching between measured temperature display and set temperature display is possible. The design is simple and easy to use.



#### Stress-free circulation connection

Connection is completed by inserting a hard tube (outer diameter of  $\phi$  10mm) only. The connector part is rotatable, so it can be turned freely and does not apply stress to the hose. Connection to a flexible hoses such as silicone hose (inner diameter of  $\phi$  9mm) is also possible by inserting the supplied



#### Safe drainage of condensation water

A drain for condensation water is provided near the connector on the upper rear side. Connect the provided drainage hose and utilize it to drain condensation water

